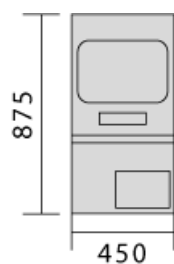


Laboratory Washer**AW47**

Model
AW47
Washing Time (Each Process)
0 to 60min. Settable
Washing Water Temp.
Room Temp. to 60deg.C

Features :

- Easy Semi-Automatic Type by turning the Switch ON after Timer Setting.
- Powerful Pressurized Water Jet from 2 directions (Top and Bottom) and from Rotating Nozzles. Detergent Washing is also available.
- Loaded with heater to make washing water hot. No boiler piping works and water heating facilities are required.
- Easy exchanging with Jet Rack (Option) for the stronger washing.

**Dimensions (mm)****Control Panel**

Specifications

Model	AW47
Performance	
Washing Type	Tow Direction, Top and Bottom, Pressurized Water Jet Type, Rotating Jet Nozzle (Exchangeable to Jet Rack)
Washing Water Temp.	Room Temp. to 60deg.C
Hot Water Supply System	Heater Heat-up (Built-in: 1kW)
Supplying Water Pressure	0.1 to 0.3Mpa
Container Stand	Table (Standard), (Racks are Option)
Water Supplying System	Electromagnetic Valve Open/Close
Water Drain Type	Natural Drainage by Level Gap
Components	
Exterior	Chrome-free Electric Galvanized Steel Plate, Chemical-proof Paint
Interior	Stainless Steel
External Dimensions	W450 x D490 x H875mm
Inner Dimensions	W420 x D450 x H570mm
Spinning Table	Φ420mm
Pump	AC100V, 50/60Hz, 200W
Door	Up/Down, Open/Close Type
Weight	Approx. 43kg
Power	Single-phase, AC100V, 50/60Hz, 150A
Accessory	
Supply Hose (with Coupler)	2m, 1pc.
Drain Hose	I.D. 25.4mm, 2m, 1pc.
Vinyl Cover, 1pc.	
Amorphous Phosphorus Detergent	1kg
Water Supply Port Unit	1 set

Optional Item (for AW47)

Jet Rack	36pcs. Settable in case of 100ml Volumetric Flask
Test Tube Rack	Approx. 450pcs. settable in case of Φ18.5mm Test Tube
Amorphous Phosphorus Detergent	8kg



Jet Rack



Test Tube Rack



Amorphous Phosphorus Detergent