

Water Purifiers


WS200/220

Water Purification Method	Distillation
Purified Water	Distilled water
Distilled Water Capacity	1.8 L/h
Collection of	Distilled water 2.2 L/min

The WS200 and 220 models are simple and low cost water purifiers especially for producing distilled water.

- The models are compact and take up little space.
- They feature a splash-prevention design to stop water from splashing and mixing with the distilled water during the boiling process. This results in stable water quality.
- The units boast a high degree of safety with the incorporation of such features as an extremely reliable device to prevent empty burning and a circuit breaker.

System Diagram

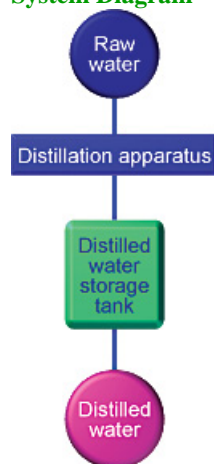

Specifications

	WS200	WS220
■ Performance		
Water purification method	Distilling method	
Purified water	Distilled water	
Distilled water capacity	Approx. 1.8 liter/Hr.	
Collection of Distilled water	Approx. 2.2 liter/min.	
■ Components		
Boiler	Stainless steel	Super hard glass
Condenser	Super hard glass	Super hard glass
Heater	Pipe heater	Set in quartz glass integument pipe
Distilled water storage tank	Polyethylene made, 20 L	
■ Attached device		
	Auto adjustment of cooling water volume, Prevention of baking heater with no load, Earth leakage breaker	
■ Standard		
Supplied water pressure range	0.5 to 5 x 100kPa (0.5 to 5kg/cm ²)	
Power source 50/60Hz	100/220V, single phase, 15/6.9A	
External dimensions(W×D×Hmm)	500×400×974	
Weight	Approx. 38 kg	
■ Accessories		
Water supply hose	2m x 1 pc. (with connecting unit)	
Drain hose	2m x 1 pc. (14mm I.D. with hose band)	
■ Consumables		
	Detergent	

■ Optional items

Remark	Description	Model No.	Product Code
○	Stand AS22	WS200/220	253176
	Supply Water Unit OWH10	WS200/220, WG35/510/710/730, WA200/500/570/710/730	253686
	Drain Trap OWI10	WS200/220, WA200	253211

Note:

The transportation cost will not be required if the order is together with main unit except the item marked "○".
In case you order the item marked "○" together with the main unit, please ask Yamato to obtain the correct transportation cost.